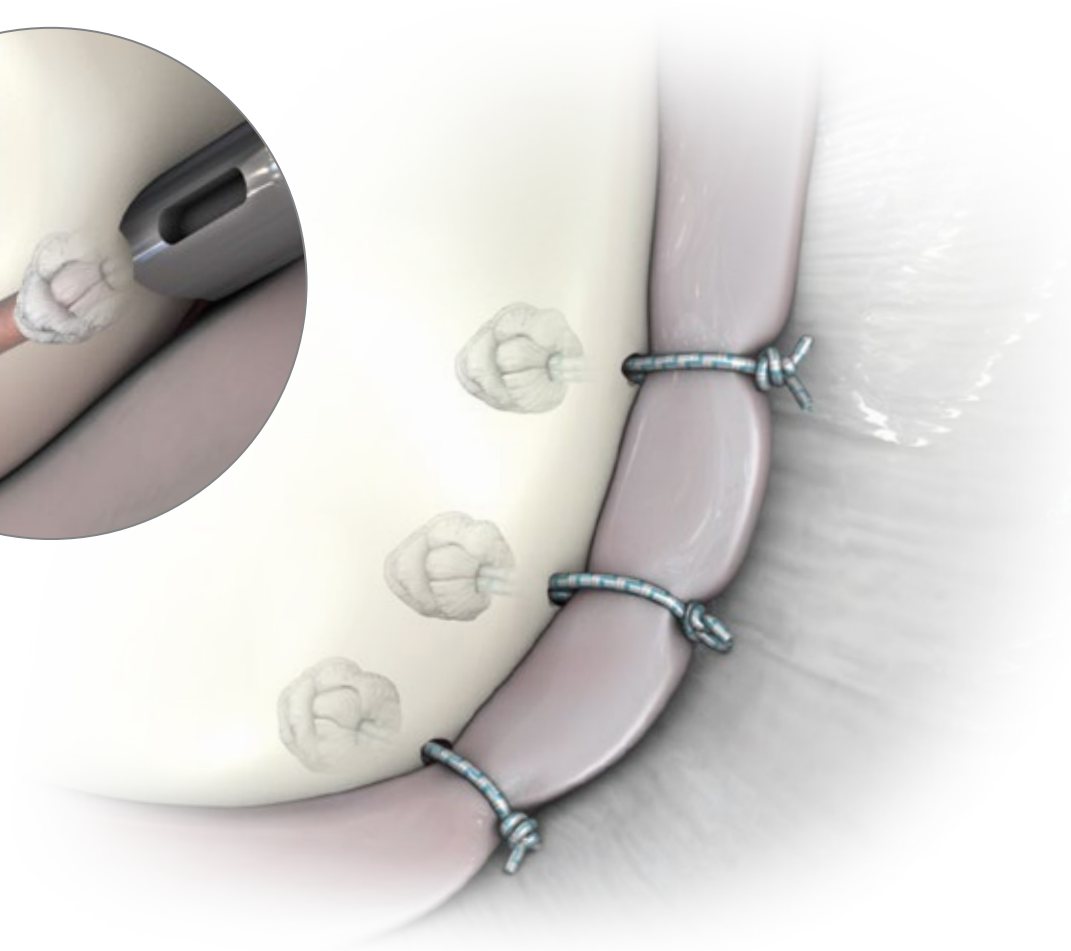
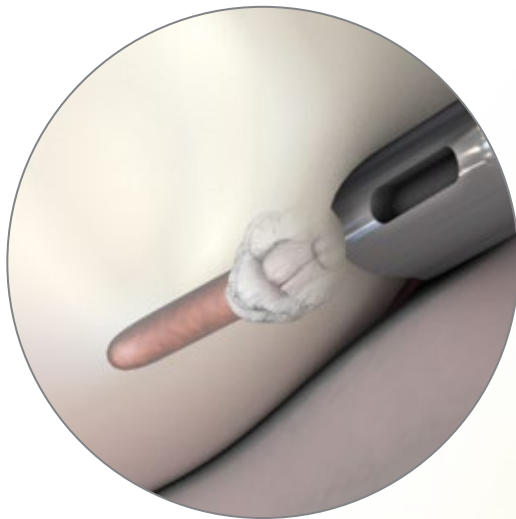


# SureLock® All-Suture Anchor System

Guide for Bankart Repair

Surgical Technique



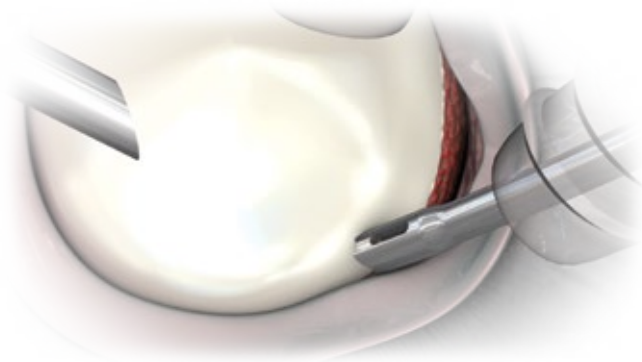


Figure 1

\*Curved and straight drill guides are available for easy access to all tear sites.

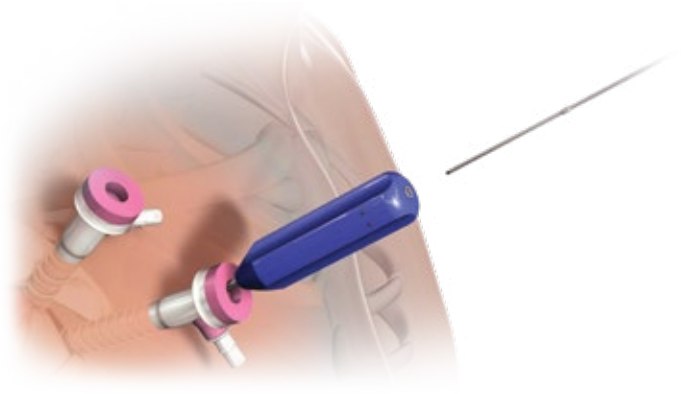


Figure 2

The Zimmer Biomet SureLock All-Suture Anchors are intended to provide soft tissue fixation to bone. The non-absorbable SureLock All-Suture Anchor is attached to an inserter and available in 1.4 mm and 2.2 mm sizes.

## Site Preparation

### Step 1

Prepare the repair site with a shaver/burr/rasp. This will help promote tendon-to-bone healing.

### Step 2

Place the SureLock obturator into the SureLock drill guide and insert through the anteroinferior portal to position the tip on the glenoid surface. Establish appropriate alignment and position of the implantation site (Figure 1).

ⓘ **Note:** The anterosuperior and anteroinferior portals are shown.

ⓘ **Note:** Curved and straight drill guides are available for easy access to all tear sites.

### Step 3

Remove the obturator and insert the SureLock drill into the proximal end of the drill guide (Figure 2).

ⓘ **Note:** The SureLock Drill Guide must be used when drilling the pilot hole.

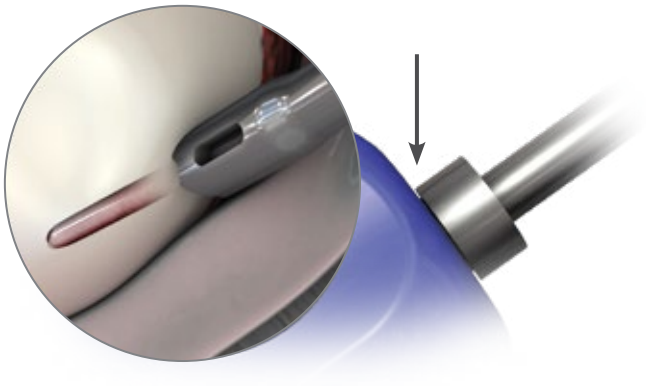


Figure 3

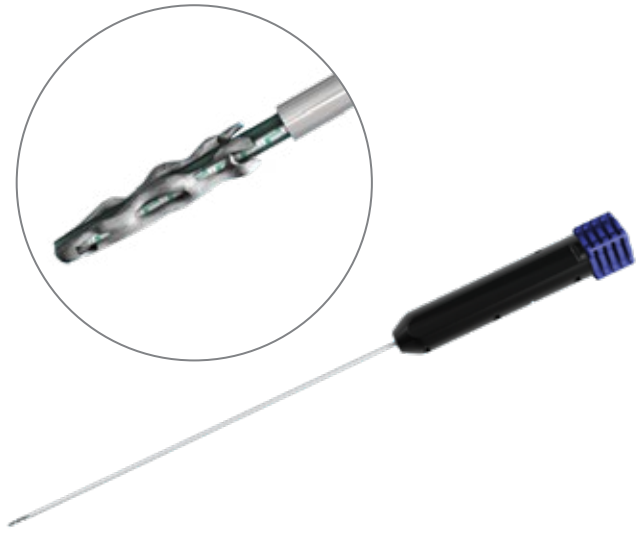


Figure 4  
SureLock All-Suture Anchor with Insertion

#### Step 4

While holding the drill guide steady on the glenoid, create a pilot hole by drilling until the stop collar on the drill contacts the proximal end of the drill guide. A positive stop will indicate proper depth has been reached (Figure 3).

#### Step 5

Remove the drill while maintaining alignment and position of the drill guide on the bone surface. This will allow for easy insertion of the SureLock All-Suture Anchor (Figure 4).



Figure 5

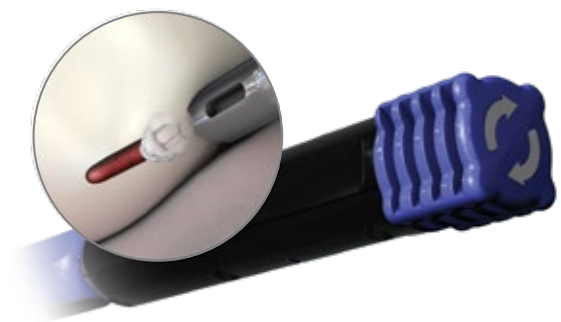


Figure 6

## Deployment

### Step 6

Insert the SureLock anchor into the proximal end of the drill guide. Lightly mallet the anchor into the pilot hole until the inserter handle comes into contact with the proximal end of the drill guide (Figure 5).

ⓘ **Note:** The SureLock Drill Guide must be used for proper deployment and fixation of the SureLock All-Suture Anchor.

ⓘ **Note:** The SureLock Anchor is fully seated when the inserter handle meets the drill guide.

### Step 7

Turn the deployment knob on the proximal end of the inserter clockwise until it stops (Figure 6).

ⓘ **Note:** Inserter-controlled delivery method eliminates manual tensioning and consistently deploys the anchor beneath the cortex.

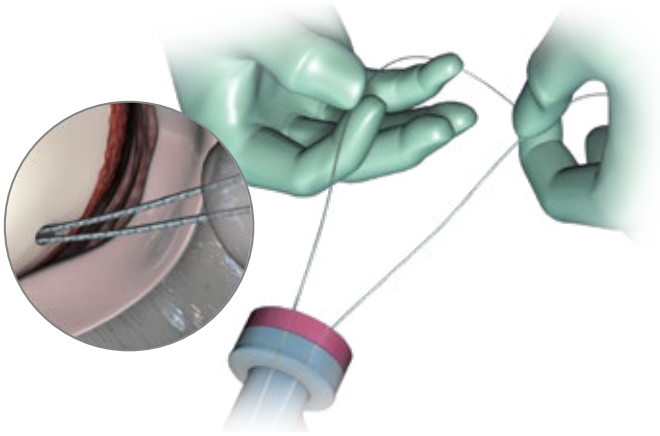


Figure 7

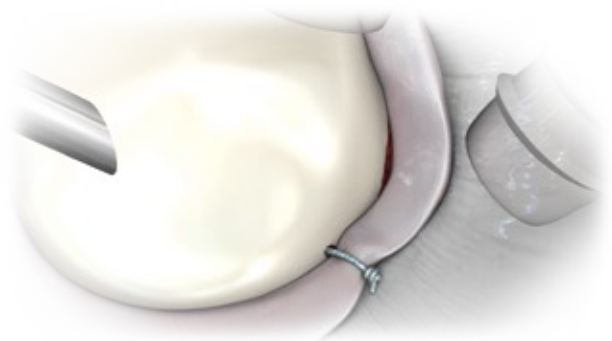
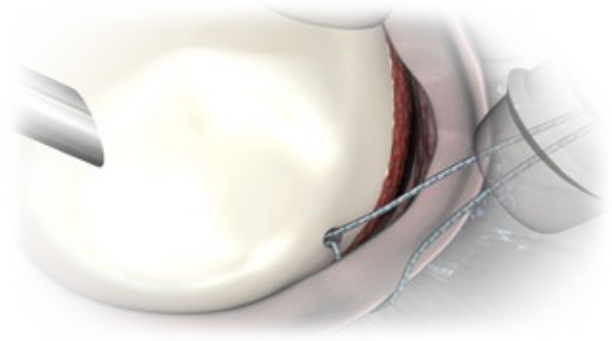


Figure 8

## Final repair

### Step 8

Pull back slowly on the inserter handle to release the inserter from the anchor. Discard the inserter handle. Remove the drill guide (Figure 7).

ⓘ **Note:** High Strength UHMWPE Suture allows for smooth, easy suture sliding and knot tying.

### Step 9

Pass suture through tissue and reattach the tissue as appropriate (Figure 8).

## Ordering Information

### SureLock All-Suture Anchor

Description	Part Number
SureLock All-Suture Anchor, 1.4 mm pre-loaded w/(1) Size 2 UHMWPE Suture, Flexible Inserter	CM-9614F
SureLock All-Suture Anchor, 2.2 mm pre-loaded w/(2) Size 2 UHMWPE Suture, Flexible Inserter	CM-9622F

### SureLock Disposable Instruments (Sterile)

Description	Part Number
Drill, for 2.2 mm SureLock All-Suture Anchor	CM-9620
Drill, with Centering Sleeve, for 2.2 mm SureLock All-Suture Anchor	CM-9620S
Drill, for 1.4 mm SureLock All-Suture Anchor	CM-9640
Drill, with Centering Sleeve, for 1.4 mm SureLock All-Suture Anchor	CM-9640S
Hard Bone Drill, 1.4 mm SureLock Anchor	CM-9640H
Hard Bone Drill with Centering Sleeve, 1.4 mm SureLock Anchor	CM-9640HS

### SureLock Instruments (Non-Sterile)

Description	Part Number
Obturator, for 2.2 mm SureLock All-Suture Anchor	CM-9621
Drill Guide (straight), for 2.2 mm SureLock All-Suture Anchor	CM-9623
Drill Guide (curved 25°), for 2.2 mm SureLock All-Suture Anchor	CM-9625
Obturator, for 1.4 mm SureLock All-Suture Anchor	CM-9641
Drill Guide (straight), for 1.4 mm SureLock All-Suture Anchor	CM-9643
Drill Guide (curved 25°), for 1.4 mm SureLock All-Suture Anchor	CM-9645

## SureLock® All-Suture Anchor

### INDICATIONS FOR USE

The Cayenne Medical, Inc. SureLock All-Suture Anchors are intended to be used for the reattachment of soft tissue to bone for the following indications:

#### Shoulder

- Capsular stabilization
  - Bankart repair
  - Anterior shoulder instability
  - SLAP lesion repairs
  - Capsular shift or capsulolabral reconstructions
- Acromioclavicular separation repairs
- Deltoid repairs
- Rotator cuff repairs
- Biceps tenodesis

#### Foot and Ankle

- Hallux valgus repairs
- Medial or lateral instability repairs/reconstructions
- Achilles tendon repairs/reconstructions
- Midfoot reconstructions
- Metatarsal ligament/tendon repairs/reconstructions
- Bunionectomy

#### Elbow

- Ulnar or radial collateral ligament reconstruction
- Lateral epicondylitis repair
- Biceps tendon reattachment

#### Hand and Wrist

- Collateral ligament repair
- Scapholunate ligament reconstruction
- Volar plate reconstruction
- Tendon transfers in phalanx

#### Hip

- Acetabular labral repair

#### Knee

- Extra-capsular repairs
  - Medial collateral ligament
  - Lateral collateral ligament
  - Posterior oblique ligament
- Patellar realignment and tendon repairs
- Iliotibial band tenodesis
- VMO advancement
- Joint capsule closure

### CONTRAINDICATIONS

- 1) Surgical procedures other than those listed in the INDICATIONS section.
- 2) Presence of infection.
- 3) Patient conditions including insufficient quantity or quality of bone or soft tissue.
- 4) Insufficient blood supply or previous infections which may hinder the healing process.
- 5) Foreign body sensitivity. If material sensitivity is suspected, testing should be completed prior to device implantation.
- 6) The use of this device may not be suitable for patients with immature bone. The physician should carefully assess the device within cartilage epiphyseal growth plates or non-ossesous tissue.
- 7) Conditions which may limit the patient's ability or willingness to follow postoperative care instructions.

This material is intended for healthcare professionals and the Zimmer Biomet sales force. Distribution to any other recipient is prohibited. For indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert or contact your local representative; visit [www.zimmerbiomet.com](http://www.zimmerbiomet.com) for additional product information. All content herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Zimmer Biomet or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet.

Check for country product clearances and reference product specific instructions for use. Zimmer Biomet does not practice medicine. This technique was developed in conjunction with health care professional. This document is intended for surgeons and is not intended for laypersons. Each surgeon should exercise his or her own independent judgment in the diagnosis and treatment of an individual patient, and this information does not purport to replace the comprehensive training surgeons have received. As with all surgical procedures, the technique used in each case will depend on the surgeon's medical judgment as the best treatment for each patient. Results will vary based on health, weight, activity and other variables. Not all patients are candidates for this product and/or procedure. Caution: Federal (USA) law restricts this device to sale by or on the order of a surgeon. Rx only.

©2018, 2022 Zimmer Biomet



**Authorized Representative**

Zimmer GmbH  
Sulzerallee 8  
8404 Winterthur  
Switzerland



**Legal Manufacturer**

Cayenne Medical Inc.  
16597 N. 92<sup>nd</sup> St  
Suite 101  
Scottsdale, AZ 85260  
USA



ZIMMER BIOMET