



Biologic Solutions

InterGro® DBM Fibers

100% DBM cortical fibers ready for
re-hydration with choice of fluid

InterGro DBM Fibers

Bone Remodeling

- 100% demineralized human cortical bone fibers offer osteoconductive and osteoinductive properties
- Each lot is verified for osteoinductivity
- Fibers provide increased osteoinductivity when compared to the osteoinductivity of particulate formulations¹



Re-hydrate with autologous blood products, saline, or antibiotic solution



Good Molding Capability²

Handling Characteristics

- User-controlled graft consistency
- Flexible and moldable upon re-hydration

Convenience

- All-inclusive bone grafting kit
- Available in 2cc, 5cc, or 10cc sizes



Ease of Preparation

Wicking Ability^{2*}



Dry

InterGro DBM Fibers



Re-hydrated

vs.



Dry

Cancellous Chips




Re-hydrated

Instructions for Use**

STEP 1

- Attach the 30cc vacuum syringe to the valve fitting on the side of the graft syringe containing DBM Fibers

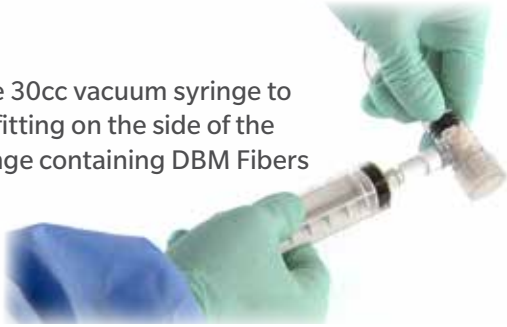

STEP 4

- Detach the dispensing syringe and discard
- Pull back and forth for 10 seconds on the plunger of the graft syringe to fully immerse the graft material with fluid



STEP 2

- Pull the vacuum syringe plunger until fully extended
- Twist plunger to engage the locking mechanism

Note: as the plunger is extended the air is being removed

- Twist off the 30cc vacuum syringe to remove



STEP 5

- Remove the cap from the end of the graft syringe
- Push the plunger gently to remove graft material

STEP 3

- Attach a dispensing syringe containing the liquid hydrating component of surgeon's choice onto the valve of the graft syringe
 - o Ensure a minimum ratio of 0.6ml fluid to 1cc DBM fibers prior to attaching dispensing syringe to graft syringe
 - o Twist the dispensing syringe in a clockwise motion to allow fluid to be dispensed into the graft syringe



Important Note: The dispensing syringe is not included in the package of InterGro DBM Fibers. Utilize an off-the-shelf luer lock syringe to add fluid.

* Each 5cc sample was hydrated with 5 ml of PRP.

** Please see package insert for complete instructions for use.

Product Availability

Product	Description	Product Number
	InterGro DBM Fibers 2cc	CDF002
	InterGro DBM Fibers 5cc	CDF005
	InterGro DBM Fibers 10cc	CDF010

To order samples contact Zimmer Biomet Spine Marketing Department.
InterGro DBM Fiber Marketing sample - CDF005S

References

1. Zimmer Biomet Internal Test Report TR-1020, Comparison of Osteoinductivity of DBM Powder to DBM Fiber
2. TR-1009, StaGraft Fiber DBM Bench Top Evaluation

INDICATIONS

InterGro DBM Fibers can be used to fill bony voids or gaps that have been surgically created, or for filling osseous defects in nonweight bearing applications.


InterGro DBM Fibers may be used with orthopedic, spinal, reconstructive, craniofacial, maxillofacial and periodontal bone grafting procedures. It can be used alone or in combination with autologous bone, or other forms of allogeneic bone in grafting procedures of non-weight bearing value. It can be used or hydrated with autologous blood, bone marrow aspirate, or autologous blood derived products such as platelet rich plasma and platelet poor plasma. It may also be hydrated with saline or antibiotic solution.

All content herein is protected by copyright, trademarks and other intellectual property rights owned by or licensed to Zimmer Biomet or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet.

This material is intended for health care professionals and the Zimmer Biomet sales force. Distribution to any other recipient is prohibited. For complete product information, including indications, contraindications, warnings, precautions, and potential adverse effects, see the package insert and www.zimmerbiomet.com.

©2018 Zimmer Biomet Spine



 **Responsible Manufacturer**
Zimmer Biomet Spine
10225 Westmoor Dr.
Westminster, CO 80021
USA

1815.2-US-en-REV0218

www.zimmerbiomet.com