




Zimmer Biomet Coonrad/Morrey Total Elbow



Interchangeability, Anterior Flange, Clinical History

The World's Number One Selling Total Elbow¹

For complete patient matching flexibility, the Coonrad/Morrey Total Elbow is unsurpassed in interchangeability and design and is the total elbow of choice for its ease of use. Backed by more than 24 years of clinical history,^{2,4} the Coonrad/Morrey Total Elbow is indicated for primary and revision intervention, post traumatic² and rheumatoid arthritis,^{3,6} and trauma reconstruction.^{4,5} The bottom line—results. That's what ease-of-operation, enhanced fixation, and simple, precise instrumentation can provide . . . that's the Coonrad/Morrey Total Elbow.

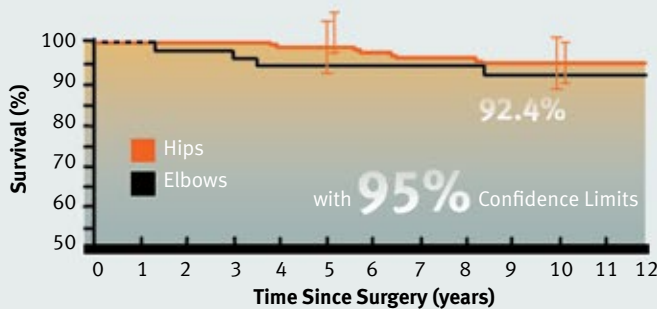

Clinical Success

The Coonrad/Morrey Total Elbow has performed successfully^{2,3} and is backed by more than 24 years of clinical history.^{2,4}





Trauma Reconstruction

The longer flange option addresses clinical situations involving distal/humeral bone loss without compromising elbow joint performance



Complete Interchangeability


Twelve Humeral Stems – 4”, 6”, and 8” lengths in extra-small, small, and regular sizes

Twelve Ulnar Stems – regular and long lengths in extra-small, small, and regular sizes

Provides comprehensive patient matching and interoperative flexibility

Two-piece locking pin simply snaps together and provides accuracy and efficiency

Easy disassembly of the locking pin enhances interchangeability in revision cases

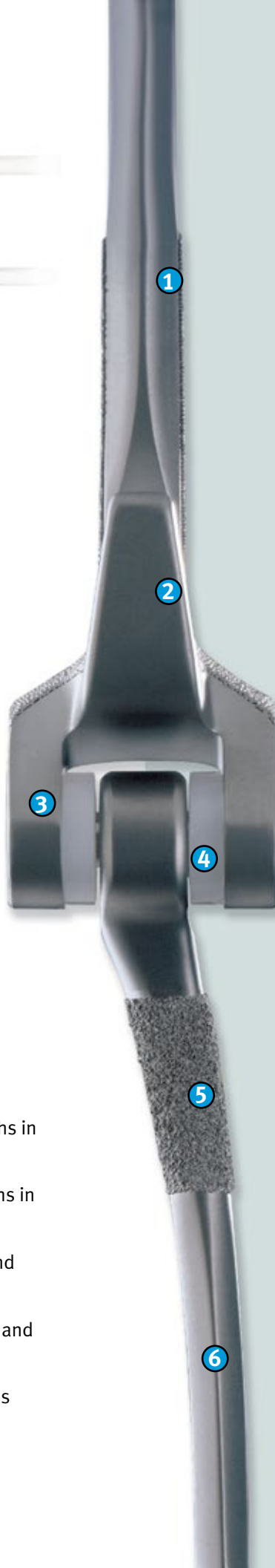

1
Triangular to match canal landscape maximizing stability by minimizing rotation

2
Proximal humeral flange permits insertion of bone graft anteriorly to enhance thickening of bone stock where maximum stress occurs, reducing posterior displacement and torsional forces and helping protect bone/cement/implant interface

3
Titanium® Ti-6Al-4V Alloy for durable, lightweight construction

4
Zimmer Poly Bushings, with ultra-high molecular weight polyethylene, prevents metal-to-metal contact

5
Plasma spray coating provides enhanced cement fixation; implant intended to be used with bone cement for both immediate and long-term fixation

6
Quadrangular to match canal landscape, to help minimize intramedullary rotation in the ulna

Implants

Humeral Implants

Prod. No.	Description
32-8105-025-04	Humeral Component, 4", Small
32-8105-025-06	Humeral Component, 6", Small
32-8105-025-08	Humeral Component, 8", Small
32-8105-026-04	Humeral Component, 4", Regular
32-8105-026-06	Humeral Component, 6", Regular
32-8105-026-08	Humeral Component, 8", Regular
32-8105-027-04	Humeral Component, 4", Extra-Small
32-8105-027-06	Humeral Component, 6", Extra-Small
32-8105-035-06	Humeral Component, Long Flange, 6", Small
32-8105-035-08	Humeral Component, Long Flange, 8", Small
32-8105-036-06	Humeral Component, Long Flange, 6", Regular
32-8105-036-08	Humeral Component, Long Flange, 8", Regular

Note: Locking pin is packaged with the humeral component.

Ulnar Implants

Prod. No.	Description
32-8105-053-01	Ulna Component, Plasma, Small, Left
32-8105-053-02	Ulna Component, Plasma, Small, Right
32-8105-063-01	Ulna Component, Plasma, Regular, Left
32-8105-063-02	Ulna Component, Plasma, Regular, Right
32-8105-073-01	Ulna Component, Plasma Long, Small, Left
32-8105-073-02	Ulna Component, Plasma Long, Small, Right
32-8105-083-01	Ulna Component, Plasma Long, Regular, Left
32-8105-083-02	Ulna Component, Plasma Long, Regular, Right
32-8105-093-01	Ulna Component, Plasma Long, Extra-Small, Left
32-8105-093-02	Ulna Component, Plasma Long, Extra-Small, Right
32-8105-043-01	Ulna Component, Plasma, 3-Inch Extra-Small, Left
32-8105-043-02	Ulna Component, Plasma, 3-Inch Extra-Small, Right

Replacement Kits* for Interchangeable System

Prod. No.	Description
32-8105-025-01	Pin/Bushing Replacement Kit Small/Regular
32-8105-027-01	Pin/Bushing Replacement Kit Extra-Small

Revision Old/New Interchange Kits*

Prod. No.	Description
32-8105-029-00	Old Regular Humeral/New Ulna Revision Kit
32-8105-030-00	Old Regular Ulna/New Humeral Revision Kit

Current Replacement/Revision Kits*

Prod. No.	Description
32-8106-000-12	Old Small Bushing Replacement/Revision Kit
32-8106-000-13	Old Regular Bushing Replacement/Revision Kit

* Consult Revision/Replacement Matrix 97-8106-004-00 for proper kit.

Provisionals

Humeral Provisionals

Prod. No.	Description
31-8105-025-04	Humeral, 4", Small
31-8105-025-06	Humeral, 6", Small
31-8105-025-08	Humeral, 8", Small
31-8105-026-04	Humeral, 4", Regular
31-8105-026-06	Humeral, 6", Regular
31-8105-026-08	Humeral, 8", Regular
31-8105-027-04	Humeral, 4", Extra-Small
31-8105-027-06	Humeral, 6", Extra-Small
31-8105-035-06	Humeral, 6", Long Flange, Small
31-8105-035-08	Humeral, 8", Long Flange, Small
31-8105-036-06	Humeral, 6", Long Flange, Regular
31-8105-036-08	Humeral, 8", Long Flange, Regular

Ulnar Provisionals

Prod. No.	Description
31-8105-053-01	Ulna, Small, Left
31-8105-053-02	Ulna, Small, Right
31-8105-063-01	Ulna, Regular, Left
31-8105-063-02	Ulna, Regular, Right
31-8105-073-01	Ulna, Long, Small, Left
31-8105-073-02	Ulna, Long, Small, Right
31-8105-083-01	Ulna, Long, Regular, Left
31-8105-083-02	Ulna, Long, Regular, Right
31-8105-093-01	Ulna, Long, Extra-Small, Left
31-8105-093-02	Ulna, Long, Extra-Small, Right
31-8105-043-01	Ulna, 3-Inch Extra-Small, Left
31-8105-043-02	Ulna, 3-Inch Extra-Small, Right

1. Data on file with Data Monitor, New York, NY 10016-5802 U.S.A.
2. Morrey BF, Adams RA, Bryan RS. Total replacement for post-traumatic arthritis of the elbow. *J Bone Joint Surg (Br)*. 1991;73-B(4):607-612.
3. Morrey BF, Adams RA. Semiconstrained arthroplasty for the treatment of rheumatoid

arthritis of the elbow. *J Bone Joint Surg*. 1992;74-A(4):479-490.

4. Morrey BF, Adams RA. Semiconstrained elbow replacement for distal humeral nonunion. *J Bone Joint Surg (Br)*. 1995;77-B:67-72.
5. O'Driscoll SW. Prosthetic elbow replacement for distal humeral fractures and nonunions. *Op Tech in Orthop*. 1994;4(1):54-57.
6. Gill DR, Morrey BF. The Coonrad/Morrey total elbow arthroplasty in patients who have rheumatoid arthritis. *J Bone Joint Surg*. 1998;89(9):1327-1335.

Instruments

Prod. No.	Description
31-8106-027-00	Humeral Guide, Extra-Small
31-8106-040-00	Hinge Pin Removal Tool
31-8106-051-00	Humeral Yolk Gauge
31-8106-052-00	Elbow Extractor
31-8106-053-01	Ulnar Rasp, Small, Left
31-8106-053-02	Ulnar Rasp, Small, Right
31-8106-063-01	Ulnar Rasp, Regular, Left
31-8106-063-02	Ulnar Rasp, Regular, Right
31-8106-065-00	Humeral Rasp, Regular
31-8106-066-00	Humeral Rasp, Small
31-8106-067-00	Starter Ulnar Rasp, Small, Left
31-8106-068-00	Starter Ulnar Rasp, Small, Right
31-8106-069-00	Pilot Ulnar Rasp
31-8106-070-00	Elbow Impactor
31-8106-076-00	Humeral Rasp, 4", Small
31-8106-077-00	Humeral Rasp, 4", Regular
31-8106-078-00	Humeral Starter Rasp
31-8106-079-00	Humeral Starter Rasp, 4"
31-8106-080-00	T-Handle
31-8106-081-00	Humeral Alignment Guide
31-8106-095-00	Humeral Cutting Guide Small/Regular
31-8106-167-00	Starter Awl, 2mm diameter
31-8106-168-00	Starter Awl, 3mm diameter
31-8106-169-00	Starter Awl, 5mm diameter
00-6601-036-00	Awl Reamer
00-8106-090-00	Elbow Instrument Case


SETS

Prod. No.	Description
31-8106-000-01	Coonrad/Morrey Instrument Set & Case
31-8106-000-02	Long Ulnar Provisional Set
31-8106-000-03	Long Flange Provisional Set
31-8106-000-04	Extra-Small Humeral Instrument Set
31-8106-000-06	Extra-Small Ulna Provisional Set
31-8106-000-07	Revision Instruments
32-8105-000-00	Standard Implant Set
32-8105-000-02	Long Ulnar Implant Set
32-8105-000-03	Long Flange Humeral Implant Set
32-8105-000-04	Extra-Small Humeral Implant Set
32-8105-000-05	Bushing Replacement Kit
32-8105-000-06	Extra-Small Ulna Implant Kit Set

This material is intended for healthcare professionals and the Zimmer Biomet sales force. Distribution to any other recipient is prohibited. For indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert or contact your local representative; visit www.zimmerbiomet.com for additional product information. All content herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Zimmer Biomet or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet.

Check for country product clearances and reference product specific instructions for use. Zimmer Biomet does not practice medicine.

The CE mark is valid only if it is also printed on the product label.

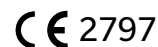

Contact your Zimmer Biomet representative or visit us at www.zimmerbiomet.com

