


# ONE Planner<sup>®</sup> Hip

Easy, Intuitive, Fast:  
Giving you a Headstart



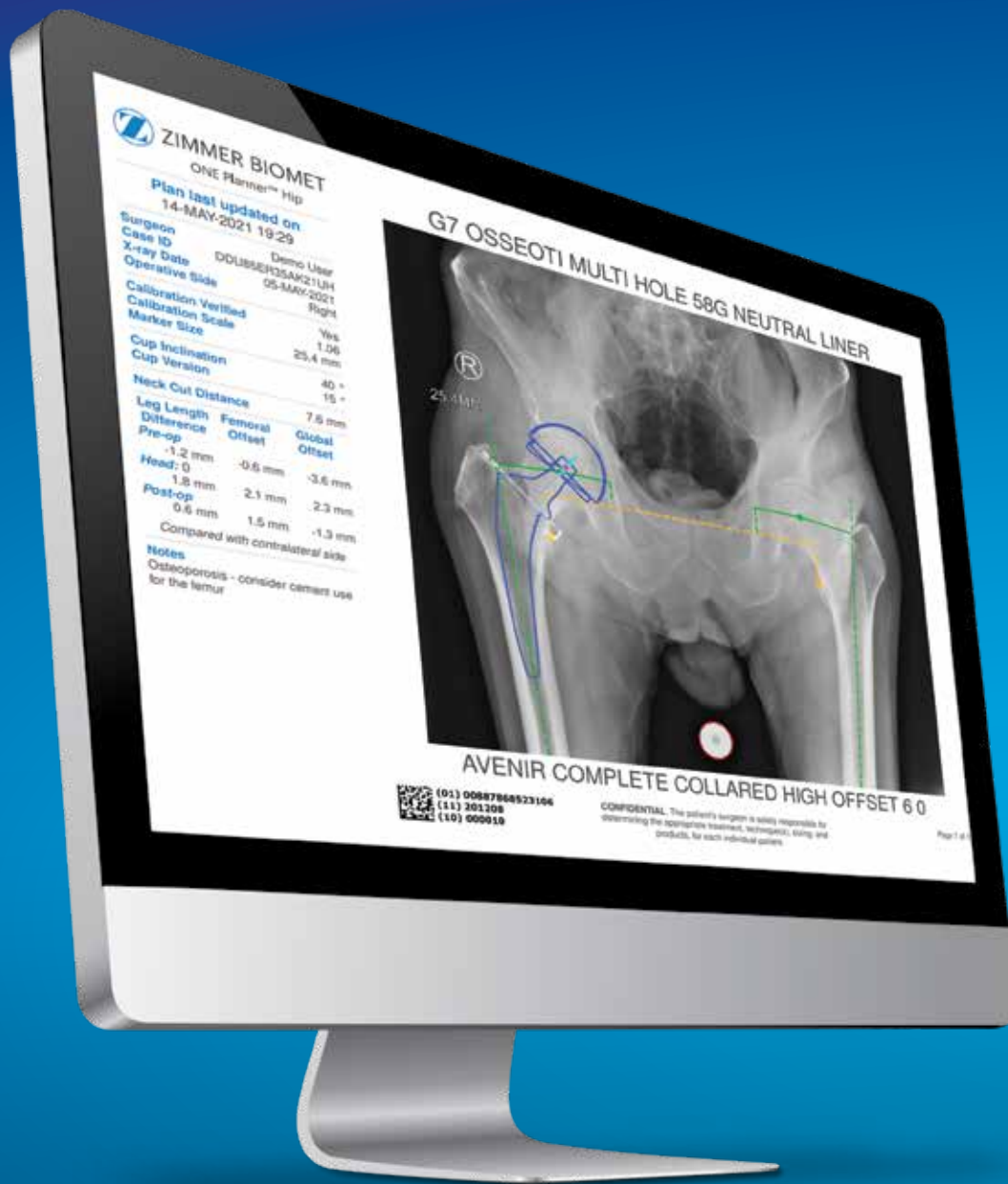

# ONE PLANNER<sup>®</sup> HIP IS AN EASY-TO-USE WEB-BASED PRE-OPERATIVE PLANNER FOR PRIMARY TOTAL HIP ARTHROPLASTY.

ONE Planner Hip requires no installation. It features a spinopelvic mobility assessment tool, if both a sitting and standing lateral X-ray are provided in addition to the AP X-ray. Auto-plan mode gives surgeons a headstart by generating a fully automated pre-operative plan, allowing surgeons to create a pre-operative plan within 5 minutes.



## EASY & INTUITIVE

- Web-based, no onsite installation needed
- No extra steps to workflow: No specific X-ray protocols are required, standard AP X-rays are suitable
- Intuitive and clean user interface with interactive adjustment of auto-plan results

## FAST AUTO-PLAN MODE

- Headstart with auto-plan mode: Pre-operative plan completed within 5 minutes
- Auto-plan mode includes image calibration, implants, offset, leg length difference and femoral neck cut

## PELVIC TILT ASSESSMENT

- Optional feature: Additional to the AP X-ray, both a sitting and standing lateral X-ray are required
- Semi-automated: Required landmarks need to be placed manually, subsequently values and angles are visualized automatically
- Categorization of patient-specific spinopelvic mobility

## TRUSTED PARTNER

- Developed in-house by the experienced Zimmer Biomet engineering team
- Compatible with Avenir Complete<sup>®</sup> Hip System, Taperloc<sup>®</sup> Complete Hip System, Echo<sup>®</sup> Bi-Metric<sup>®</sup> and the G7<sup>®</sup> Acetabular System
- ONE Planner Hip is suitable for use with the ROSA<sup>®</sup> Hip system

## HOW TO ACCESS

Access via existing Zimmer Biomet case management systems: SMS for U.S., DCMS for outside the U.S.




This material is intended for health care professionals. Distribution to any other recipient is prohibited.

For indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert or contact your local representative; visit [www.zimmerbiomet.com](http://www.zimmerbiomet.com) for additional product information.

All content herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Zimmer Biomet or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet. ©2021, 2023 Zimmer Biomet

Check for country product clearances and reference product specific instructions for use.



 **Legal Manufacturer**  
Zimmer CAS  
75 Queen Street Suite 3300  
Montreal (Quebec) H3C 2N6  
Canada