



www.ldrmedical.com

LDR, LDR Spine, LDR Medical, Avenue, BF+, BF+(Ph), Bi-Pack, C-Plate, EasySpine, Laminotome, L90, MC+, Mobi, Mobi-C, Mobi-L, Mobidisc, ROI, ROI-A, ROI-C, ROI-MC+, ROI-T, SpineTune and VerteBRIDGE are trademarks or registered trademarks of LDR Holding Corporation or its affiliates in France, the United States and other countries.

United States
13785 Research Boulevard
Suite 200
Austin, Texas 78750
512.344.3333

France
Hôtel de bureau 1
4, rue Gustave Eiffel
10 430 Rosières Prés Troyes, FRANCE
+33 (0)3 25 82 32 63

China
Unit 06, Level 16, Building A,
Beijing Global Trade Center #36
North Third Ring Road
East, Dongcheng District,
Beijing, China, 100013
+86 10 58256655

Brazil
Av. Pereira Barreto, 1395-19º Floor
Room 192 to 196
Torre Sul-Bairro Paraiso
Santo André / São Paulo-Brazil
CEP 09190-610
+55 11 4332 7755



Avenue[®] L

LATERAL LUMBAR CAGE with



Avenue® L

LATERAL LUMBAR CAGE

The Avenue® L lateral lumbar cage is designed for strength, versatility and stability. It represents the pinnacle of lateral lumbar cage evolution incorporating zero-profile, intradiscal, and integrated fixation. The enhanced in-line, self-guided VerteBRIDGE® Plating Technology facilitates simplified cage insertion all through a direct, minimally invasive approach.

Features

- Comprehensive offering of footprints to meet different patient anatomy and provide for optimal endplate contact
- Bevelled nose to ease insertion
- I-beam design to increase the rigidity of the cage and support the graft during insertion into the intervertebral space
- Intra-operative adjustment of the Cage Holder to ensure optimal positioning of the cage prior to insertion of VerteBRIDGE plating
- Guided in-line plate delivery
- Sterile packaging for assured product quality and absolute traceability


Dimensions

Product Range and Sizing †

Widths	17mm and 22mm
Lordotic Angles	0° and 6°
Heights	8mm, 10mm, 12mm, 14mm and 16mm*
Lengths	40mm, 45mm, 50mm, 55mm and 60mm
Plate Lengths	Short, Medium and Long

† Product availability may vary by country and market.

* Special Order

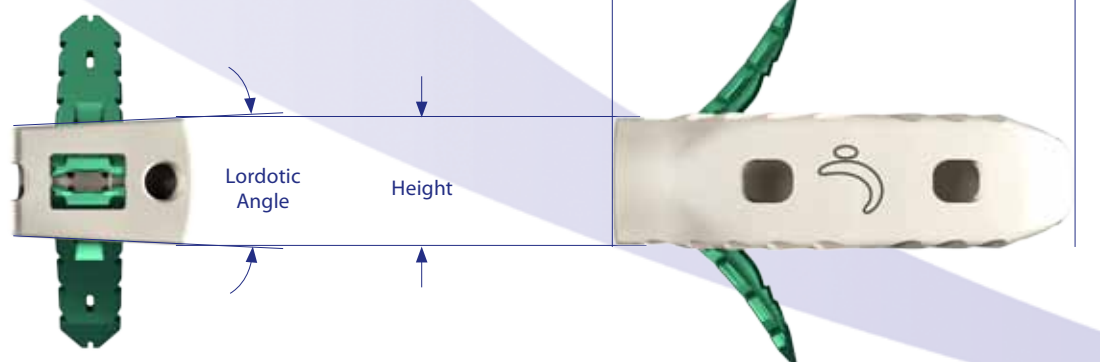

Featuring VerteBRIDGE® Plating Technology

- **Zero-profile** design with no hardware protruding lateral of the vertebral bodies
- **Self-locking** plates designed for initial and long-term stability
- **Self-guided, curved** plates ease insertion
- Inserter provides **protection of anatomic structures** during plate placement
- **Innovative plate delivery system** optimized for the direct lateral approach



Indications for use in the U.S.:

The Avenue® L Interbody Fusion System is indicated for intervertebral body fusion of the lumbar spine, from L2 to S1, in skeletally mature patients who have had six months of non-operative treatment. The device is intended for use at either one level or two contiguous levels for the treatment of degenerative disc disease (DDD) with up to Grade I spondylolisthesis. DDD is defined as back pain of discogenic origin with degeneration of the disc confirmed by history and radiographic studies. The device system is designed for use with or without integrated fixation and must be used in conjunction with posterior supplemental fixation (e.g. pedicle screws). The device system is intended to be used with autograft to facilitate fusion.