



ZIMMER BIOMET

Your progress. Our promise.®



Cervical Solutions

**Trellss<sup>™</sup>-C**  
Porous Ti Interbody System

A new foundation for growth

# TrellOss™ -C

## Porous Ti Interbody System

A 3D printed titanium interbody platform featuring a scaffold structure with 70% porosity and a 7 micron roughened surface topography to foster a cellular relevant environment for adhesion and bone ingrowth.<sup>1</sup>


### TrellOss-C Implant

- Rigid teeth help to resist implant migration
- Central window for graft packing and containment
- Implants are sterile-packed for reduced risk of contamination and hospital reprocessing costs
- Zero-profile inserter for access and visualization of disc space
- Removable depth stop for inserter/trials to accommodate surgeon preference

### TrellOss-C Sizes

HEIGHTS	FOOTPRINT	LORDOSIS
5 mm–12 mm	12 mm x 14 mm	6°   0°
5 mm–12 mm	14 mm x 16 mm	6°   0°

# A NEW FOUNDATION FOR GROWTH

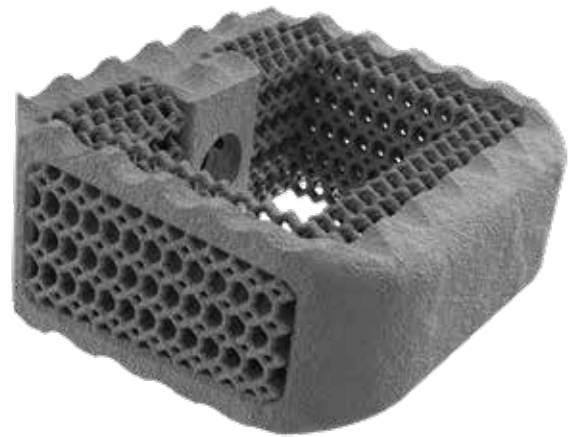

## **Porosity**

Open architecture with 70% porosity including varying pore sizes of 300, 500, and 700 microns that mimic cancellous bone allowing for a conducive environment for cellular activity<sup>1,5,6,7</sup>

## **Texture**

7 micron surface texturing creates an environment for potential cellular adhesion<sup>2,3,4</sup>

## **Structure**

Scaffolding structure provides additional surface area<sup>2,3</sup>

## References

1. McGilvray KC, Easley J, Seim HB, et al. Bony ingrowth potential of 3D-printed porous titanium alloy: a direct comparison of interbody cage materials in an in vivo ovine lumbar fusion model. *Spine J* 2018;18(7):1250-1260.
2. Olivares-Navarrete R, Hyzy SL, Slosar PJ et al. Implant materials generate different peri-implant inflammatory factors: poly-ether-ether-ketone promotes fibrosis and microtextured titanium promotes osteogenic factors. *Spine* 2015;40(6):399-404.
3. Olivares-Navarrete R, Hyzy SL, Gittens RA, et al. Rough titanium alloys regulate osteoblast production of angiogenic factors. *Spine J* 2013;13(11):1563-70.
4. Rao PJ, Pelletier MH, Walsh WR, et al. Spine Interbody Implants: Material Selection and Modification, Functionalization and Bioactivation of Surfaces to Improve Osseointegration. *Orthop Surg* 2014;6:81-89.
5. Ponader S, von Wilmsowsky C, Widenmayer M, et al. In vivo performance of selective electron beam-melted ti-6al-4v structures. *J Biomed Mater Res A* 2010;92A:56-62.
6. Li JP, Habibovic P, et al.: Bone ingrowth in porous titanium implants produced by 3D fiber deposition. *Biomaterials* 2007;28:2810.
7. Karageorgiou V, Kaplan D. Porosity of 3D biomaterial scaffolds and osteogenesis. *Biomaterials* 2005;26(27):5474-91.



Manufactured by: Nexxt Spine, LLC  
14425 Bergen Blvd, Suite B  
Noblesville, IN 46060



**ZIMMER BIOMET**  
Your progress. Our promise.®

800.447.3625/zimmerbiomet.com

©2019 Zimmer Biomet Spine, Inc. All rights reserved.

All content herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Zimmer Biomet Spine, Inc. or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet Spine. This material is intended for health care professionals, the Zimmer Biomet Spine sales force and authorized representatives. Distribution to any other recipient is prohibited.

For product information, including indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert and [www.zimmerbiomet.com](http://www.zimmerbiomet.com). Check for country product clearances and reference product specific instructions for use.  
2634.1-GLBL-en-REV0619